

## Product Datasheet

**Product code:** **394**

**Product type:** UltraCut Scissors

**Made in:** Italy


**Description:** The 394 are UltraCut scissors with fine blades, specifically designed for miniature works such as cutting suture fibers. The short cutting edges combined with the extra long handles boost the leverage, making cutting easily and effortlessly. The sharp blades guarantee a precise and clean cuts. The micro-serrated blade allows non-slip cutting. The special Inox material combined with the polish finishig is the best solution for cleaning. The blades are hardened to 60HRC hardness.

**Tip description:** straight, fine, sharp, micro-serrated

**Coating/finishing:** Mirror Polished

**Body material:** Stainless Steel (Hardened AISI 440)

**Picture:**

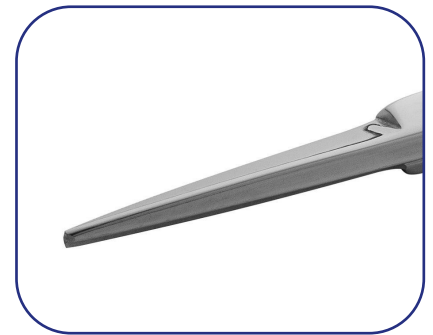

**Product dimensions:**

|        |                       |
|--------|-----------------------|
| OAL    | 112.0 mm – 4.409 inch |
| Width  | 46.0 mm – 1.811 inch  |
| Height | 4.5 mm – 0.177 inch   |

**Tip dimensions:**

|               |                      |
|---------------|----------------------|
| Tip width     | 0.70 mm – 0.028 inch |
| Tip thickness | 0.70 mm – 0.028 inch |